

## Demag\_GUI Usage and Tips (Version 3.2)

### Launching

The Gui may be launched either with the command line by going to the directory containing demag\_gui.py and running it with:

```
python ./demag_gui.py
```


If you have QuickMagIC already running you can launch demag\_gui by just clicking on the demag\_gui button in the main window shown below:




### Interpretation of Specimen Data

#### Adding Interpretations:

You can analyze specimen component data by adding fits with the add fit button. Additionally you can select the fit you would like to edit or view by using the drop down box under the add fit button. Once you have selected a fit the shape of the selected fit's data points will turn to a diamond shape to distinguish them from the other data points.

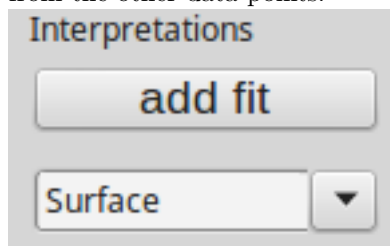

Once you have selected the desired fit you can edit its bounds using the drop down boxes under the bounds header

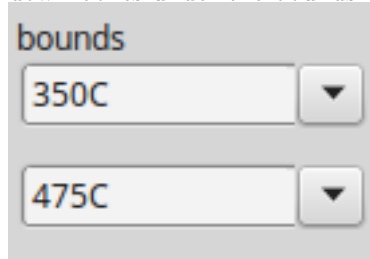

Alternatively you can double click the list of measurement steps on the left to pick out the bounds of your interpretation. The included steps in the currently selected interpretation are shown in blue on this list and measurement steps marked bad are shown in yellow. (This method of selecting bounds is recommended for mac users who have in the past experienced problems with the drop down boxes)

i	Step	Tr	Dec	Inc	M
7	T	450.0	179.1	-0.0	4.83e-0
8	T	480.0	169.9	9.7	4.54e-0
9	T	500.0	169.2	12.5	4.21e-0
10	T	510.0	165.2	14.6	3.97e-0
11	T	520.0	159.2	15.7	3.66e-0
12	T	530.0	164.3	16.7	3.26e-0
13	T	540.0	165.3	16.7	3.02e-0
14	T	545.0	160.8	18.6	2.48e-0
15	T	550.0	159.7	15.8	2.11e-0
16	T	555.0	170.6	17.2	1.64e-0
17	T	560.0	165.5	18.2	1.30e-0
18	T	565.0	162.0	18.5	8.99e-0
19	T	568.0	161.3	19.8	6.30e-0
20	T	571.0	153.7	25.0	3.39e-0
21	T	571.0	162.0	20.8	3.27e-0
22	T	574.0	175.0	11.5	1.86e-0
23	T	577.0	121.0	67.1	1.14e-0
24	T	577.0	130.0	57.1	9.11e-0

You may notice that the fit will be given a generic name such as *Fit 1* you can change the name of the fit from default by typing into the drop down box containing fits then pressing enter. You can anchor your interpretation or preform a plane fit using the drop down box under specimen mean type (default: line).

specimen mean type

line ▼

Coordinate Systems available as well as orientation of Zijderveld are available on the left.

specimen:

Z35.1a ▼

previous next

---

coordinate system:

geographic ▼

---

Zijderveld plot:

X=NRM dec ▼

Specimen data for a fit can be seen in the upper center of the GUI in a large box labeled specimen mean statistics.

specimen mean statistics

dec	inc	n	mad	dang	a95
177.9	-44.8	5	3.2	17.8	

### Flagging Bad Measurement Data

You can set acceptance criteria to a pmag\_criteria table by using Analysis/“Acceptance Criteria”/“Change Acceptance Criteria”. If any measurement steps are bad you can flag them as such by right clicking on the list of measurement steps to the left of the GUI. If you flag a step bad that you would later like to restore you can simply right click on it again and it will be flagged as good again.

## Saving Specimen Interpretations

Once you have picked out your interpretations you can save the session data in two different ways a .redo file or pmag tables. In addition you may save image files of the plots.

**The .Redo File:** You can use Analysis/“Save current interpretations to a redo file” to create this file type or you can just hit the save button next to add fit, this method is recommended as it prevents accidental pressing of the clear all interpretations button. **Note:** this file type does **NOT** load previous interpretations on start up you must go to the menu option Analysis/“Import previous interpretations from a redo file” to restore your previous session.


**The Pmag Tables:** By going to the menu File/“Save MagIC pmag tables” you can export your interpretations made in Demag GUI to the MagIC pmag tables which can then be used by other MagIC programs or uploaded to the MagIC database. You can export any or all of the 3 coordinate systems upon selecting this option and you may choose to save pmag\_samples, pmag\_sites, and pmag\_results tables in addition to the pmag\_specimens table that is output. If you choose to output additional information you will be prompted by a pop up window for additional information. **Note:** this save format loads on start up of the GUI immediately restoring your session, also selection of this option will overwrite your demag\_gui.redo file in the working directory.

**Images of Plots:** Select the menu option File/“Save plot”/“Save all plots” to save all plots alternatively you can save any of the plots individually. Some examples can be seen below:


# Zijderveld plot

X-axis rotated to NRM (177); NRM=3.93e-07 Am<sup>2</sup>

Z35.1a

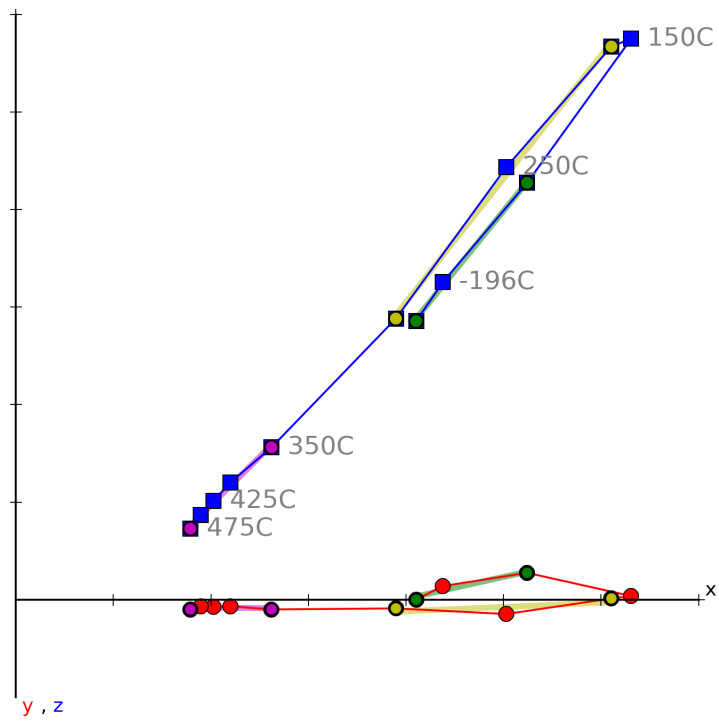


specimen: Z35.1a

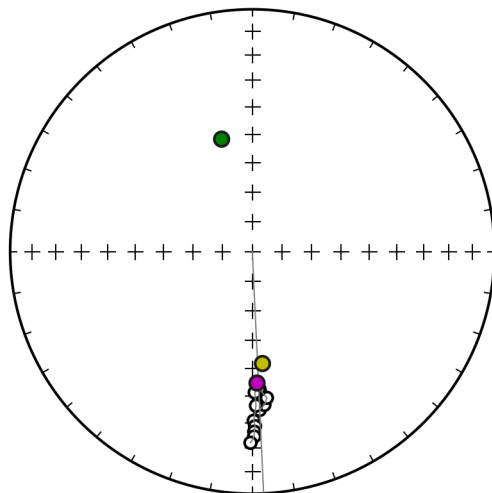


M/M0

Z35.1a

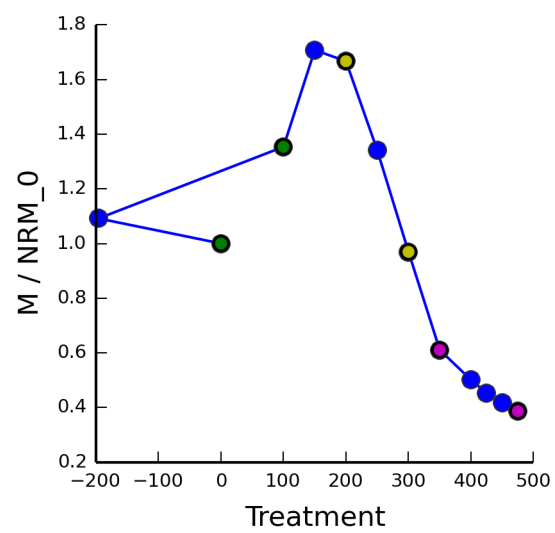


site: Z35

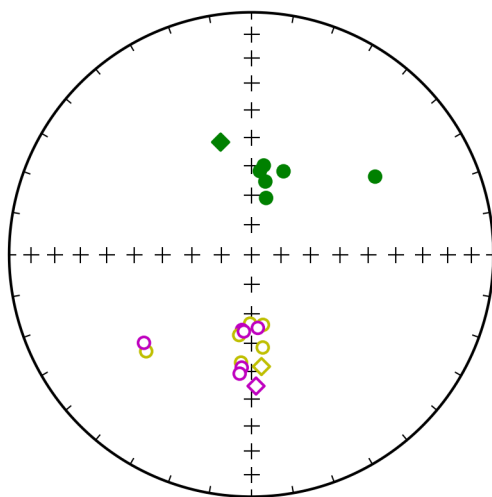


## Deleting Specimen Interpretations

If you would like to delete a single interpretation select the one you wish to delete from the interpretation drop down box then click delete. Alternatively if you wish to clear all interpretations you may go to the menu option Analysis/“Clear all current interpretations”.

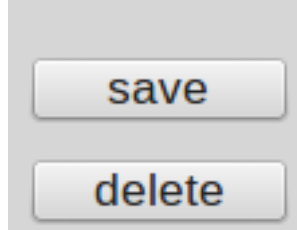


## Higher Level Plots and Interpretation

The set of drop down boxes to the right of the interpretation data is there to determine what level you want to analyze in the higher level analysis options include: site, sample, location, and study. The drop down below this selects which of the available sites, samples, location, or studies you want to display.

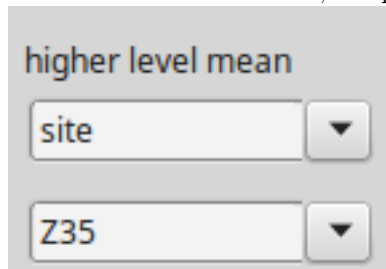

You can then select how to group your data by using the drop down menu under the show header. You can select what kind of mean to take using the first drop down under the mean header. Which interpretations to use for the means can be selected under the second drop down menu. You can then use the remove/replace button to remove or replace the set of points belonging to the current specimen in the higher order mean.

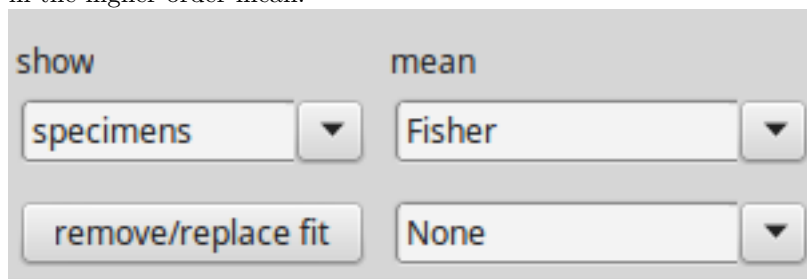

You can view the higher order stats results in the bottom left of the GUI.

**Fisher statistics:**

dec: 183.9

inc: -54.0

alpha95: 48.1

K: 1.0

R: 6.8791

n\_lines: 21

n\_planes: 0